



Protocube Manual

Functions and Definitions

Fluorescence Microscope






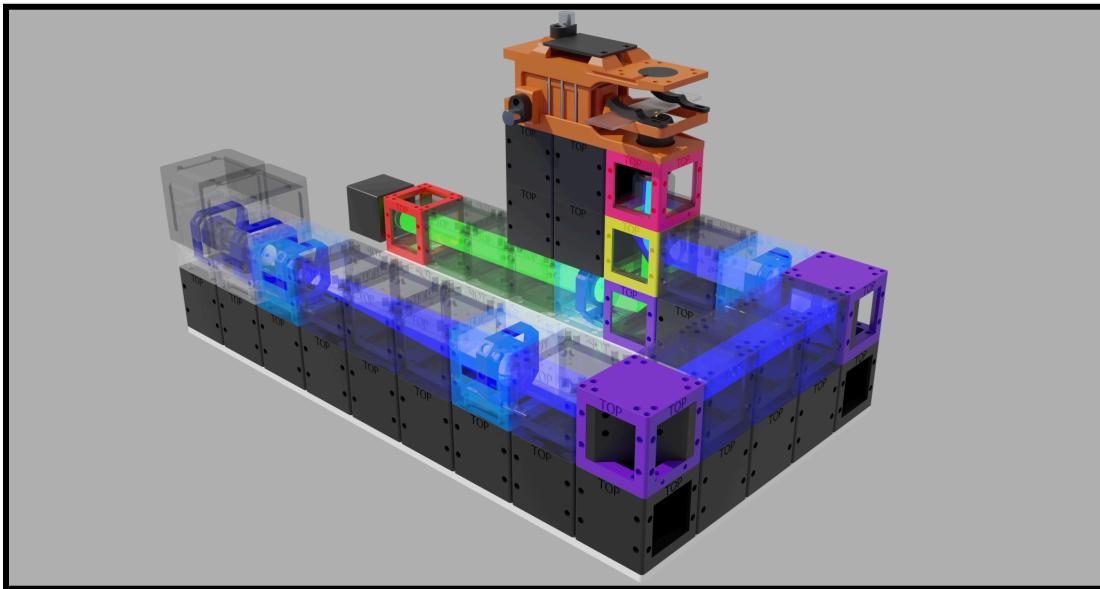


Transmitted Light Microscope



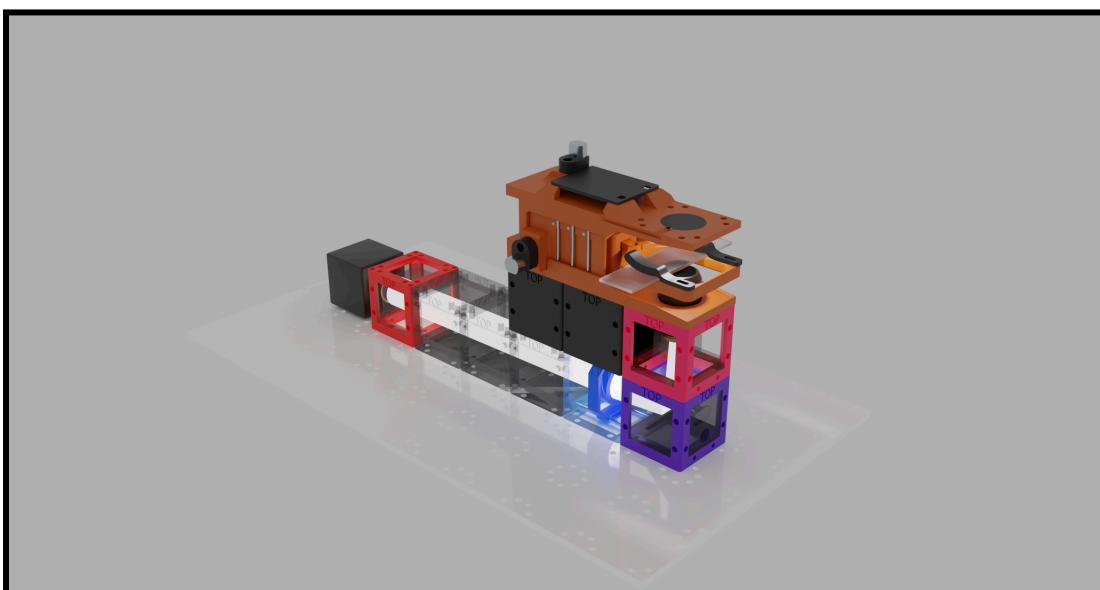


Protocube Manual

Cubes!

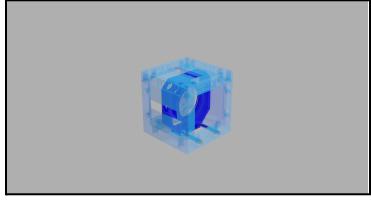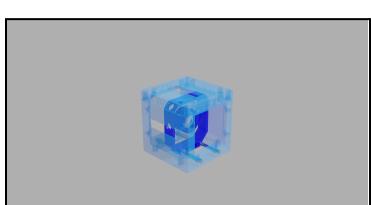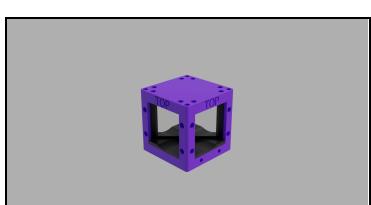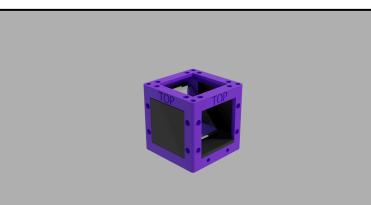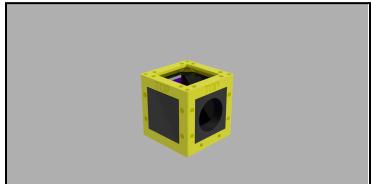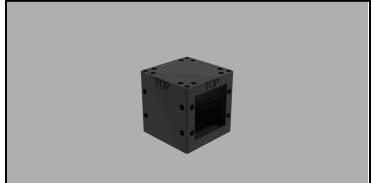Image	Name	Description
	One Inch Lense Cube	A Light Blue Cube designed to hold a standard 1" optical lens, the lens is fixed in place by an M4 set screw and a press fit cap for extra security. The lens is aligned positionally by two M6 set screws and a manual rail adjustment.
	Half Inch Lense Cube	A Light Blue Cube designed to hold a standard 1/2" optical lens, the lens is fixed in place by an M4 set screw and a press fit cap for extra security. The lens is aligned positionally by two M6 set screws and a manual rail adjustment.
	Horizontal Mirror Cube	A Purple cube designed to fit a rectangular mirror. The mirror's distance and alignment are determined by three brass thread inserts heat pressed in, along with an M3 Screw contacting magnets on the mirror plate.
	Vertical Mirror Cube	A Purple cube designed to fit a rectangular mirror. The mirror's distance and alignment are determined by three brass thread inserts heat pressed in, along with an M3 Screw contacting magnets on the mirror plate.
	Dichroic Cube	A Yellow cube designed to hold a standard 1" Dichroic. It is sized for SM1 to hold standard 1" filters via press fit.
	Basic Cube	A Black cube used for support and light occlusion. It prevents unnecessary noise and signal interruption and provides an easy solution to a safely enclosed light path.

Protocube Manual

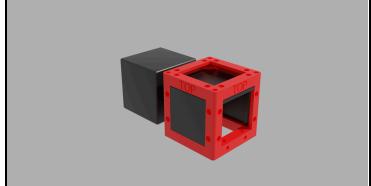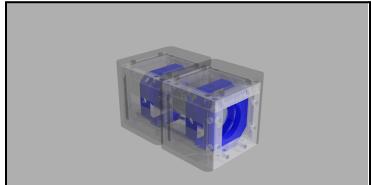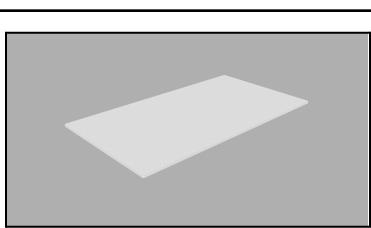Image	Name	Description
	Camera Cube	A Red cube designed to fit a standard size SM1 via a press fit and keep the light path enclosed.
	Objective Cube	A Pink cube designed to hold a standard microscope objective. It is placed at the appropriate distance on the z-axis for use with a stage.
	Light Source Cube	A White cube that holds the light source and is adjustable in the same way as the lenses, but with two mounts to align and tilt the light source. It also includes covers to be enclosed once aligned.
	XYZ Stage!	An XYZ-Stage designed to use consumer electronics from Pimoroni and put them together into a functional, position control stage. Currently the encoders are unused, and other attachments for the stage are available.
	Baseplate	A laser-cut Base Plate that holds magnets arranged to help align and snap together the cubes.

Protocube Manual

Assembly Instructions

What you need:

- At least 24x Black Basic Cubes and 1x fully open Basic Cube or more
- At least 2x Horizontal and 1x Vertical Mirror or more
- At least 1x Light Source Assembly and a fitting Dichroic Cube or more
- At least 1 x Half Inch Cube f30 and 3x One Inch Cubes f200 or more
- At least 2x Targets, one of each type or more
- At least 1x Objective Cube
- At least 1 x Camera Cube
- At least 1x Stage



SAFE USE GUIDANCE - GENERAL

A Class 3R laser is low powered. It normally would not harm eyes during a momentary exposure of less than ¼ second. This is within the aversion response, where a person turns away and/or blinks to avoid bright light.

Do not deliberately look or stare into the laser beam. Laser protective eyewear is normally not necessary. A Class 3R laser is not a skin or materials burn hazard.

However, a Class 3R laser can be a distraction, glare or flash blindness hazard for pilots and drivers. **NEVER aim any laser towards an aircraft or vehicle that is in motion.** This is unsafe and is illegal -- you could be arrested and jailed.

ONLY ALLOW USE BY RESPONSIBLE PERSONS!!!


This is **NOT a toy**. Children can safely use Class 3R lasers only with continuous adult supervision.

Laser Source	Power (With ND Filter)	Power (Without ND Filter)
Red (640nm)	13.9mW	109.2mW
Blue (405nm)	7.9mW	49.7mW
Green (532nm)	4.3mW	29.5mW

Reminders!

- Never touch any lenses or mirrors directly with your hands or fingers!
- Never put any part of your body in the beam path, it is dangerous!

Protocube Manual

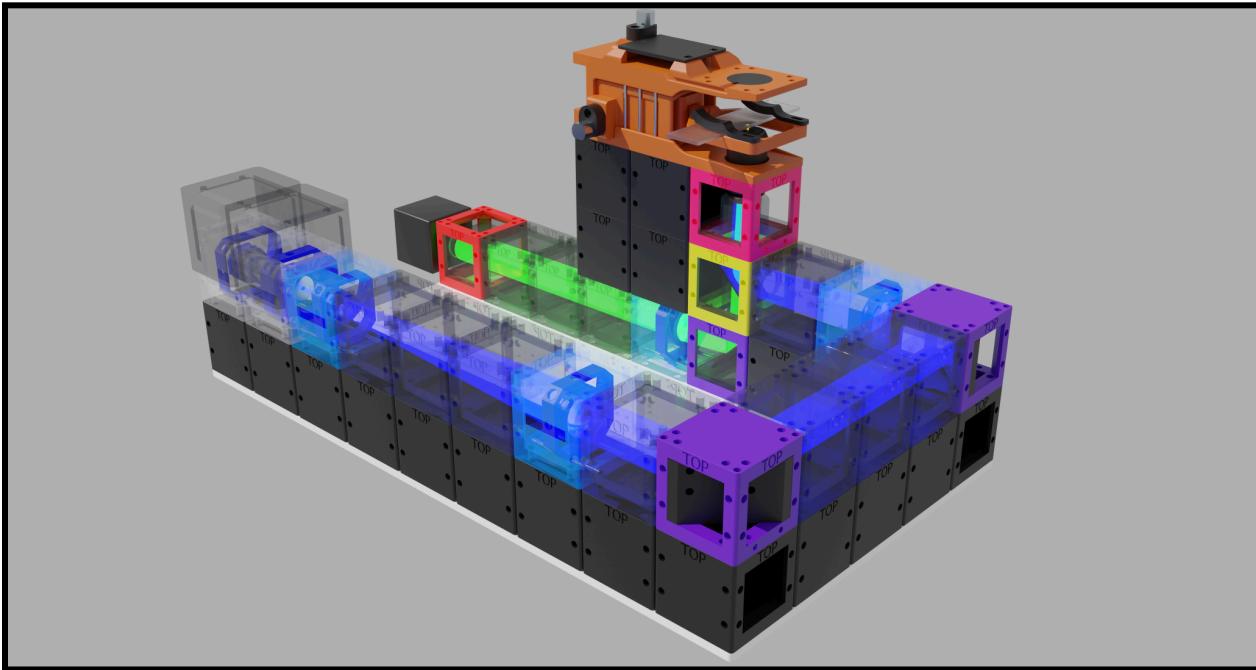

Excitation Beam Path:

- **Laser Alignment**

1


- Start by placing a laser source on top of two of the black basic cubes, placing them on top of the magnetic base plate. Remove the laser source covers and loosen the four set screws by unscrewing them a little bit.

2

- Place two targets in front of the laser like in the image below, making sure to also place them on top of the magnetic base plate at least three cubes apart from each other, and one cube away from the laser.

3

- Set up a flat target such as paper on the other end of the targets, to see how the laser needs to be adjusted.

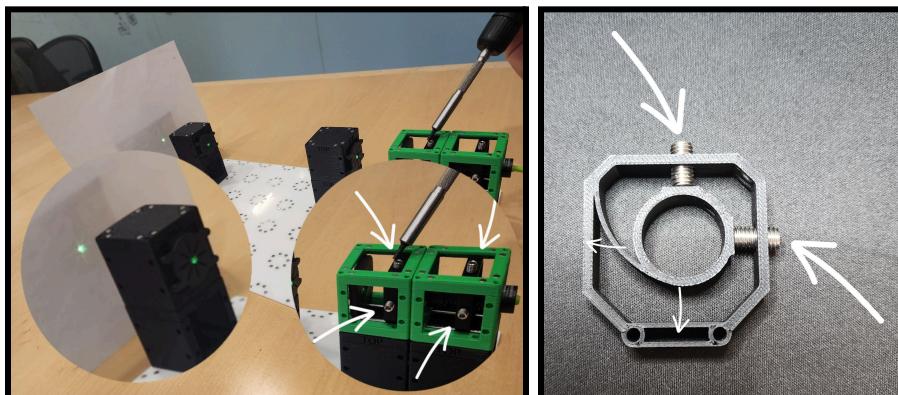


Protocube Manual

4

- The set screws push the inner mount horizontally and vertically to allow for precise control of movement.
- Using this, you can make small adjustments, centering and tilting the laser to shine through the center of both targets, swapping them around to double check that it is aligned correctly like above.

- **Recommendations**

- Tap the laser source cube on the top with your screwdriver after every adjustment to make sure the internal mount is seated correctly.
- Slightly tilt the cubes, as the magnets don't set exactly into the same spot and can be adjusted by wiggling the cube a tiny bit.

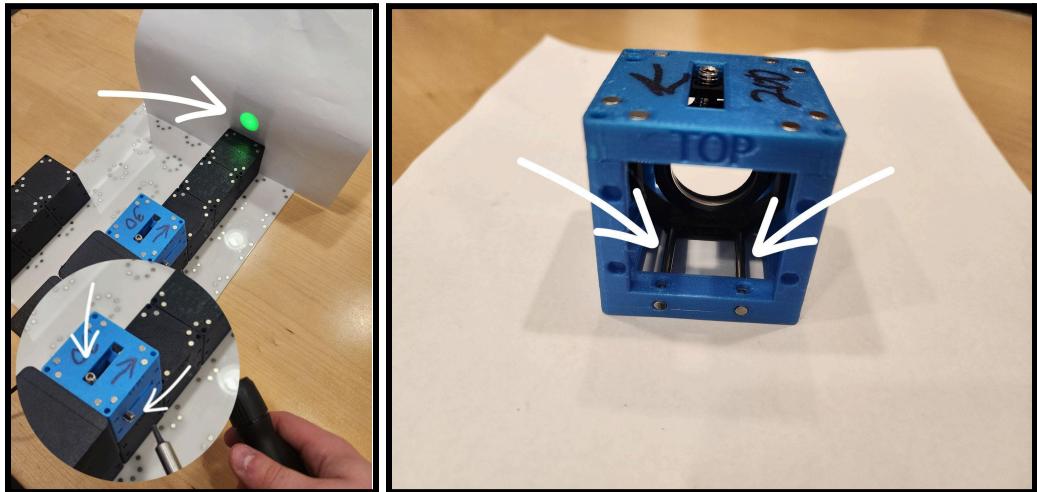

- **Beam Collimation**

1

- Start by placing an f30 lens directly in the laser path, loosening the set screws by unscrewing them a little bit.
- Slide the lens back on its rails to be as close to the laser as possible, this is easier done from the bottom of the cube, while picked up.
- Use the set screws to center it as close as possible to the beam path like in the image above. Placing a cube behind the paper with light behind it to act as a target to center to.

2

- Next place an f200 lens about 200mm away and center it by using the same method as before with the f30 lens.


Protocube Manual

3

- Slide the f200 lens back and forth along its rails in the cube until it is focusing the beam path to be aligned and parallel, or collimated, and the beam coming out doesn't shrink or grow when moving a flat target such as paper farther away.

4

- You may have to move the f30 lens back and forth as well to get this as best as you can.

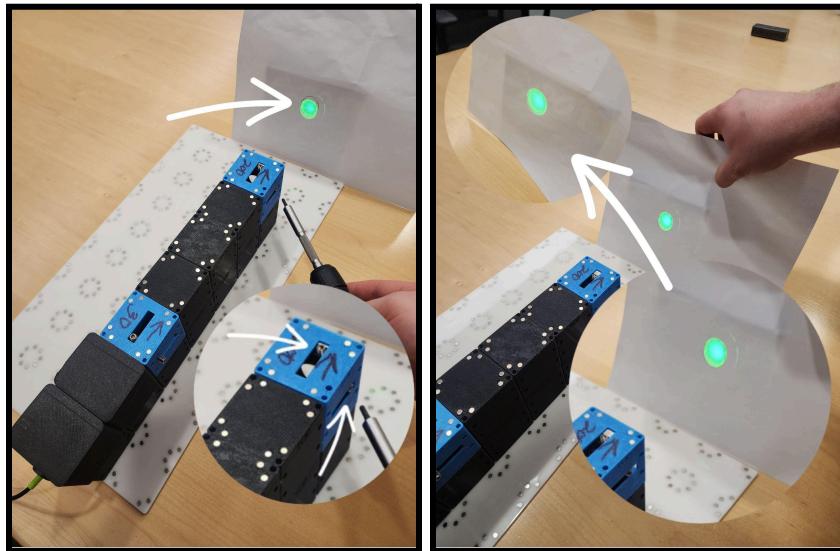


● Mirror Alignment

1

- If necessary, direct the beam with the horizontal mirrors, making sure it stays in the same plane by adjusting the screws on the back to tilt the mirror and adjust its distance like in the images below.

2

- I recommend using the same paper method used in the beam collimation steps to make sure that the beam path stays straight like in the images below.

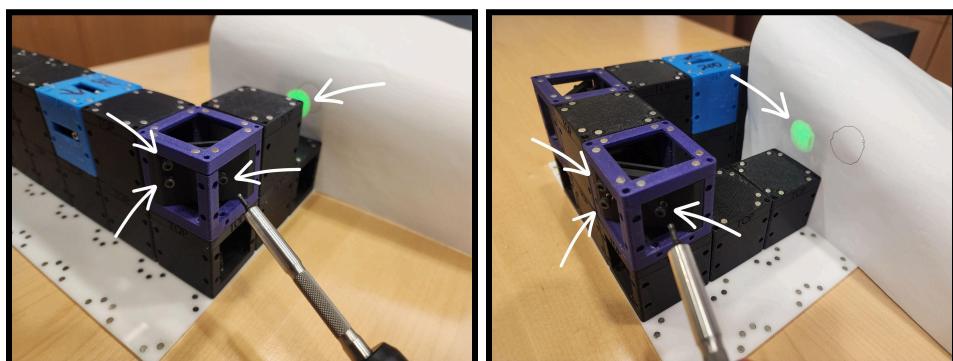

Protocube Manual

- **Objective Focusing**

1


- Direct the beam into an f200 lens, making sure that it is centered with the set screws.

2


- Focus the beam into the objective like a cone using the f200 lens, sliding it back and forth so that the beam path comes out on top of the infinity corrected objective parallel.

3

- Direct the beam into a correct set of filters and dichroic mirror in the dichroic cube, and up into the objective.

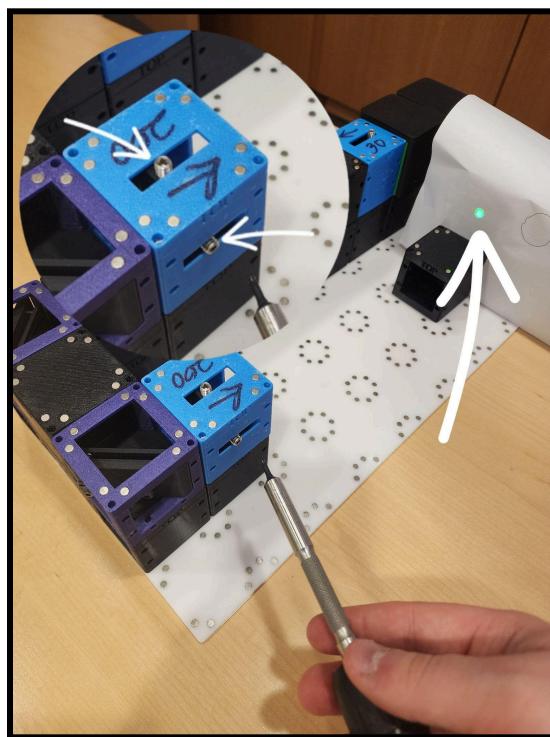

Protocube Manual

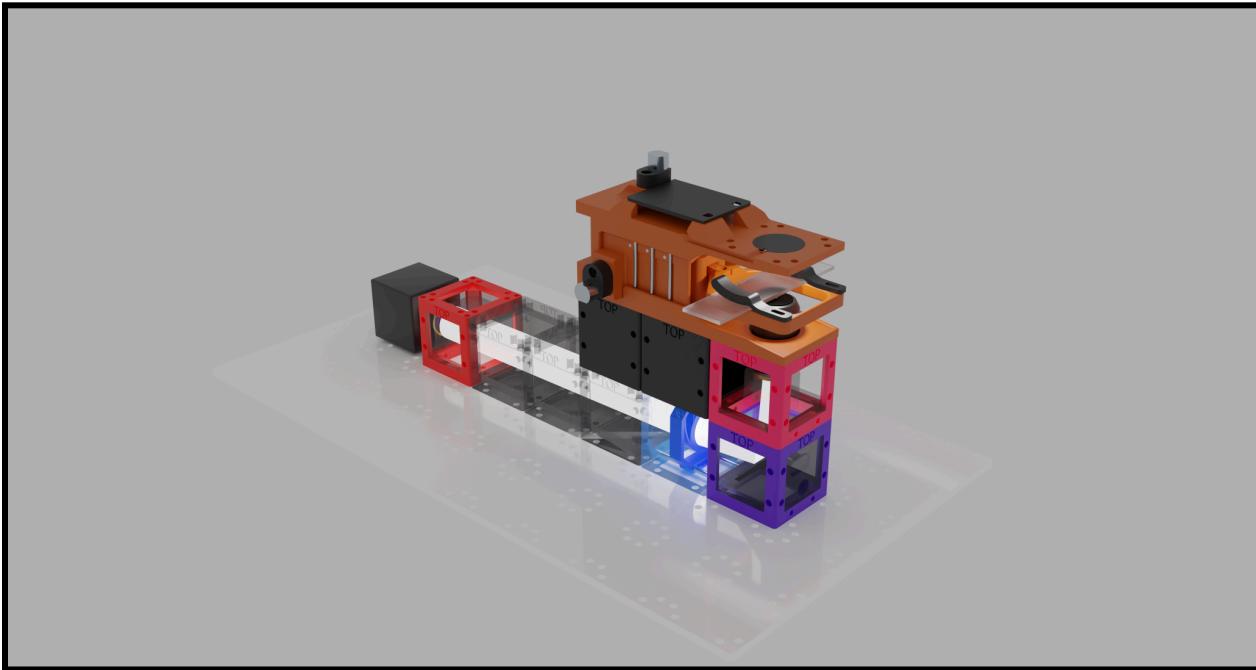


Emission Beam Path:

- 1**
 - Place the XYZ-stage on top of the objective.
- 2**
 - Let the beam pass through the filter cube and onto a vertical mirror, adjusting the mirror via its screws to keep the beam parallel and direct it forward.
- 3**
 - Place a f200 lens in front of the beam, making sure it's centered using the set screws.
- 4**
 - Place the camera so the camera sensor is about 200mm from the f200 lens. Start focusing on the stage

Protocube Manual

Gamepad

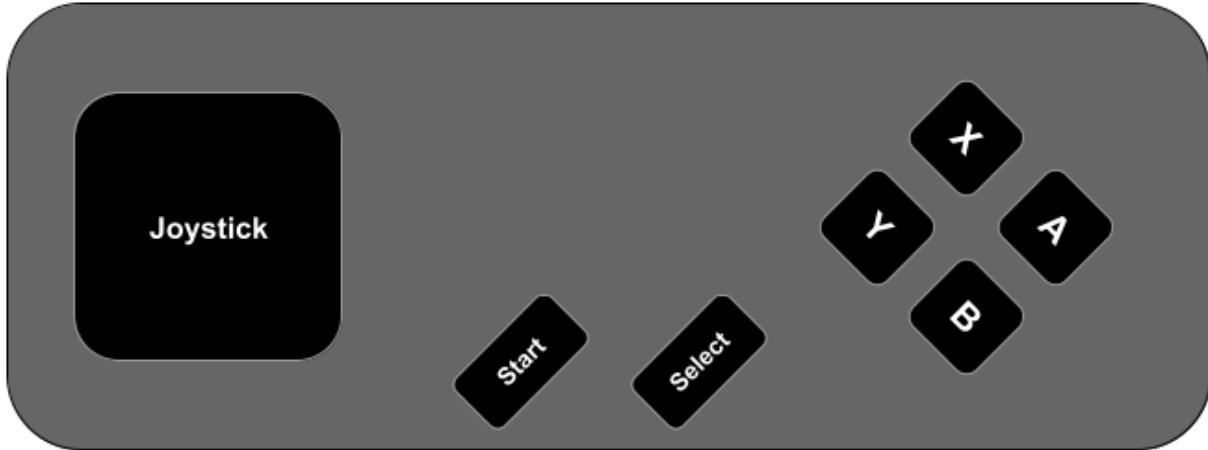

Joystick

Used to move the X, Y, and Z Axis

(It does not contribute to speed control, push it all the way in one direction at a time to move the stage)

Start

Switch between X/Y, and Z mode on the Joystick

(Hold for one second, Z Mode only uses the vertical axis of the joystick)

Select

Reset any false positive limits, click if you're having issues.

(Hold for one second, PAY ATTENTION TO THE MECHANISM, IT WILL TEAR ITSELF APART If You're not careful)

Y

Motor Speed 27.5%

X

Motor Speed 62.5%

A

Motor Speed 100%

B

Toggle LED On/Off

(Hold each button for a second to confirm the speed or setting is set correctly)